

COMPRESSION STOCKINGS

A Guide for Patients

+oapl.™



INTRODUCTION

You have been fitted with a pair of graduated compression stockings. Compression stockings provide support and aid in venous (blood) return. The stockings are graduated to move fluid only toward the heart. The stockings compensate for the inefficiency of your veins or lymphatic system in moving fluid and blood back to your heart.

Compression stockings come in various compression ratings and styles to cope with different conditions and severity. Your Orthotist will have measured you and fitted you with the correct size and style according to your specific requirements.

APPLICATION

Stockings should be fitted first thing in the morning and removed before going to bed at night.

There are several ways to put stockings on and your orthotist will show you the best way for you.

Opposite is the most common method:

Consulting Orthotist

Clinic Attended

Date Attended

FITTING




Place your hand into the top of the stocking and run your fingers down to the heel. Pinch the heel with your forefinger and thumb.

With your other hand pull the leg section of the stocking inside out whilst holding the heel in your fingers.

With the heel under your foot draw the stocking foot onto your foot until it reaches the front of your ankle. Ensure that the heel is in the correct place.

Gather the leg section of the stocking up and fold back over the toes and gradually pull up the leg evenly. Avoid twisting the material, pulling it at the edges, or stretching the stocking along its length.

Note: If you have overstretched the stocking push it back down to the ankle, then pull it up again a little at a time, not stretching it too far.

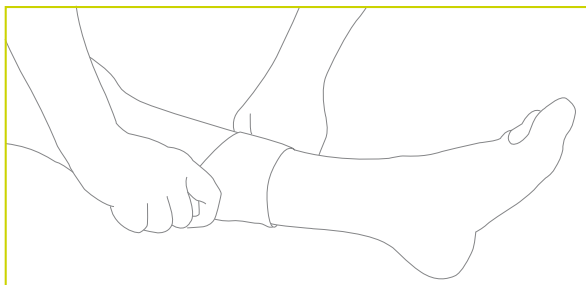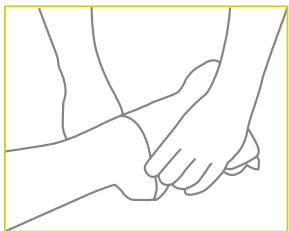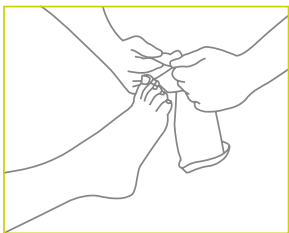

FITTING

Check that the top band of the stocking is smooth and flat.

Under no circumstances is it to be rolled or folded over .

Smooth out all wrinkles by running your palm over the stocking.

Note: Below knee stockings should sit about 3cm below the crease in the back of your knee.



This process may be made easier by;

Placing a plastic bag over your foot before pulling the stocking on (this is only possible for open toe garments).

Wearing rubber gloves (used for washing dishes) to help you grip the stocking.

- For moist skin use talcum powder to help the stocking glide over the skin.

- Ask for assistance from a relative or a friend!

Also available are silky slippers (instead of plastic bags) and special frames to aid in putting stockings on. Your orthotist will suggest these options should you require them.

REMOVAL

The easiest method of removing a stocking is to do the reverse of the way that you put it on.

Many people have trouble moving the stocking past their heel, as this is the tightest section, but pushing your thumb down the back of your heel and easing it over your heel will make it easier.

GENERAL CARE

Stockings are best washed by hand in warm water with mild soap (eg wool wash) and rinsed in clear cold water. Do not use bleach. Roll stocking into a towel and gently squeeze excess water out. Dry in the shade or indoors. Avoid sunlight or direct heat. The stockings must not be ironed or tumble dried. Do not hang stocking from toe or top as it may stretch, instead lie it flat on a dry towel.

Stockings are best washed each evening before bed and should dry by morning. However it is recommended to have a spare stocking to wear on alternate days.

Although stockings are quite strong be careful of sharp finger nails and rings! To avoid damaging the stocking use rubber gloves when putting the stocking on.

Stockings should always be worn with shoes to protect the stocking and prevent you from slipping.

To increase durability, socks can be worn over the top of the stockings.

Ensure that footwear is in good condition on the inside since a defect in the shoe lining may increase friction and quickly lead to damage of your stocking.

Please DO NOT repair the stocking yourself. If a defect occurs contact your orthotist.

Over time stockings begin to stretch and will no longer provide adequate support. At this time they will probably need replacing.

We usually recommend that patients be remeasured for subsequent / replacement pairs of stockings as your legs may have less swelling or sizing may have changed.

GENERAL CARE

Compression stockings should only be worn after consultation with your doctor. If any of the following signs appear remove the stocking immediately and contact your orthotist:

- Swelling of the toes or above where the stocking finishes.
- “Blue” toes
- Change in sensation in your foot (ie tingling, numbness)
- Pain
- Laddering of the fabric (in new stocking)

CONTACT DETAILS

Fitzroy Clinic

93 - 97 Webb Street
Fitzroy VIC 3065
T: 03 9419 2499
F: 03 9416 3534
mainfacility@fitzroy.oapl.com.au

Box Hill Clinic

Epworth Eastern
Suite 5A, Level 2, Arnold Street
Box Hill VIC 3128
T: 03 9897 1374
F: 03 9783 6944

Bundoora Clinic

Northpark Private Hospital
Greenhills Road
Bundoora VIC 3083
T: 03 9467 1222

Canberra Clinic

Mobility Matters
35 Townsville Street
Fyshwick ACT 2609
T: 02 6280 7244
F: 02 6239 1281
oapl@mobilitymatters.com.au

Clayton Clinic

281 Clayton Road
Clayton VIC 3168
T: 03 9562 9422
F: 03 9562 9605
clayton@oapl.com.au

Footscray Clinic

55 Pickett Street
Footscray VIC 3011
T: 03 9687 4589
F: 03 9687 3961
footscray@oapl.com.au

Frankston Clinic

346 Nepean Highway
Frankston VIC 3199
T: 03 9783 3866
F: 03 9783 6944
frankston@oapl.com.au

Malvern Clinic

1489 Malvern Road,
Malvern VIC 3144
T: 03 9824 5874
F: 03 9824 8205
malvernclinic@oapl.com.au

Richmond Clinic

Epworth Centre
Level 7, Suite 5, 32 Erin Street
Richmond VIC 3121
T: 03 9421 6226
F: 03 9426 4321

Ringwood Clinic

86 Mt Dandenong Road
Ringwood VIC 3135
T: 03 9879 4299
F: 03 9870 9411
ringwood@oapl.com.au

Bendigo Clinic

Anne Caudle Centre
Barnard Street
Bendigo VIC 3550
T: 03 5454 8752
F: 03 5454 8756
oapl@bendigohealth.org.au

Main Facility and Administration

93-97 Webb Street

Fitzroy Victoria 3065

T: 03 9419 2499

F: 03 9416 3543

mainfacility@oapl.com.au

ACN 006 520 314

ABN 30 006 520 314

www.oapl.com.au

Version 4.09

OAPT0019

